

ULTRASONIC WAVE TOOTHBRUSH



3 BRUSHING MODES

Solution for different oral problem



Normal Clean



Soft



Massage

- Toothbrushes can be used by children and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Cleaning and user maintenance shall not be made by children without supervision.
- Children shall not play with the appliance.
- The supply cord cannot be replaced. If the cord is damaged the appliance should be scrapped.
- The appliance is only to be used with the power supply unit provided with the appliance.

RLT206 ULTRASONIC WAVE TOOTHBRUSH IMPORTANT SAFEGUARDS

The RLT206 has been carefully designed to offer you and your family a new and unique brushing experience, that is both safe and effective. However, When using the electrical products, especially when children are present, basic safety precautions should always be followed, including the following:

READ ALL INSTRUCTIONS BEFORE USING DANGER TO REDUCE THE RISK OF ELECTROCUTION:

- Always unplug this product immediately after using.
- Do not use while bathing.
- Do not place or store product where it can fall or be pulled onto a tub or sink.
- Do not place in or drop into water or other liquid.
- Do not reach for a product that has fallen into water unplug immediately.
- Specific instructions for the proper method of cord storage, total product storage, and the like when the product is not in use; and for cord while in use, such like.

WARNING

To reduce the risk of fire, electric shock, or injury:

- Do not leave Appliance when plugged in. Unplug from outlet when not in use and before servicing.
- Do not allow to be used as a toy. Close attention is necessary when used by or near children.
- Use only as described in this manual. Use only manufacturer's recommended attachments.
- Do not use with damaged cord or plug. If Appliance is not working as it should, has been dropped, damaged, left outdoors, or dropped into water, return it to a service center.

- Do not pull or carry by cord,use cord as a handle,close a door on cord,or pull cord around sharp edges or corners.
- Do not run Appliance over cord.Keep cord away from heated surfaces.
- Do not unplug by pulling on cord.To unplug,grasp the plug,not the cord.
- Do not handle plug or Appliance with wet hands.
- Turn off all controls before unplugging.
- Never block the air openings of the product or place it on a soft surface,such as a bed or couch,where the air openings may be blocked.Keep the air openings couch,where the air openings may be blocked.Keep the air openings free of lint,hair,and the like.
- Never use while sleeping or drowsy.
- Never drop or insert any object into any opening or hose.
- Do not use outdoors or operate where aerosol(spray)products are being used or where oxygen is being administered.
- Connect this product to a properly grounded outlet only.See grounding instructions.
- Unplug this product before filling.Fill(reservoir)with water only unless otherwise specified by manufacturer.
- Do not overfill(or specify filling instructions).

SAVE THESE INSTRUCTIONS

Household USE Only


PLEASE READ CAREFULLY BEFORE USE:

- This product has no user serviceable parts.
- Store the charging unit onto a dry and safe place when it is not in use.
- Untwisted the power cord before use.
- In case of the product failing to function,you should not dismantle,exchange or manipulate any part of it,otherwise there is risk of an electric shock.You should take it to the seller.
- If the supply cord is damaged,it must be replaced buy the manufacturer in order to avoid a hazard.You should contact the seller.
- If the product is dropped,The brush head should be replaced before the next use even if no damage is visible.
- This appliance has a polarized plug(one blade is wider than the other).As a safety feature,this plug will fit in a polarized outlet only one way.If the plug does not fit fully in the outlet,reverse the plug.If it still does not fit,contact a qualified electrician. Do not attempt to defeat this safety feature.

Declaration

The toothbrush is conformity to the applicable provisions of applicable US Federal requirements,including 21 CFR,Chapter 1,Subchapter J,Radiological Health.

LIST OF COMPONENTS

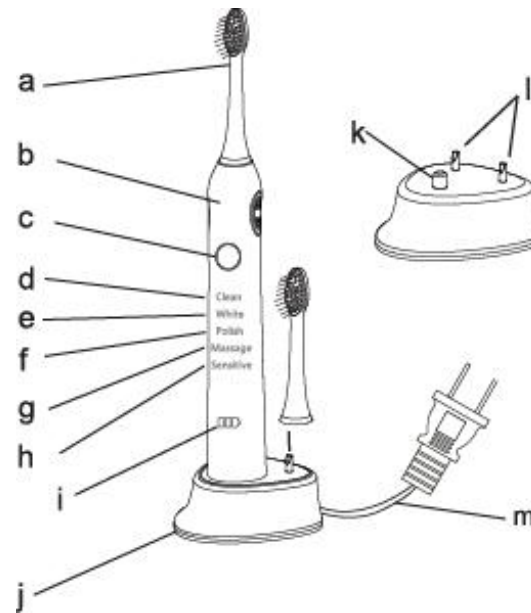

- a. Brush head
- b. Toothbrush Handle
- c. On/Off & Brushing Mode Button
- d. Clean Working Mode Indicator
- e. Whiten Working Mode Indicator
- f. Polish Working Mode Indicator
- g. Massage Working Mode Indicator

- h. Sensitive Working Mode Indicator
- i. Power & charge indicator
- j. Charger Unit
- k. Handle Stander
- l. Brush head Stander
- m. Cable

CONNECTING AND CHARGING

- RLT206 has a waterproof handle. It is electrically safe and is designed for use in the bathroom. But it's not meant to be used while underneath the water.
- A full charge takes at least 24 hours and enables approximately 6 weeks brushing. The constant blue LED ("Ch") indicates charging is in progress. Once fully recharged, the light turns off automatically.
- Pull plug out of outlet for saving energy, or leave it and keep the brush on charging base, there is no risk of overcharging of battery.

USING THE RLT206

Attach a brush head to the handle. You can use any kind of toothpaste. To avoid splashing, guide the brush-head to your teeth before switching the appliance on.

BRUSHING TECHNIQUE:

- Guide the brush head slowly from tooth to tooth.
- Hold the brush head in place for a few seconds before moving on to the next tooth.
- Brush the gums as well as the teeth. Start with the outside of the teeth, then the insides ending with the chewing surfaces.
- Do not scrub or press hard, simply let the brush do all the work.

BRUSHING MODES

The RLT206 offers five effective brushing modes for varying oral care needs. Features a auto shut-off after two minutes of brushing and a dentist-recommended auto-alert when brushing the same area for more than 30 seconds.

Clean -with up to 31,000 strokes per minute ,exceptional mouth cleaning

White -with up to 41,000 strokes per minute ,removes microscopic plaque and bacteria left behind by other toothbrushes

- | | |
|------------------|--|
| Polish | -With up to 48,000 strokes per minute ,2X more plaque will be effectively scrubbed away so fast |
| Massage | -with up to 31,000 strokes per minute,strong to weak,weak to strong in a circle , gentle stimulation of gums |
| Sensitive | -with weak 31,000 strokes per minute ,gentle cleaning for sensitive areast |

TWO MINUTE AUTO-TIMER:


The RLT206 is designed to reach the proven and best brushing timing.After 2 minutes of brushing,the device stops in 2 mins automatically.

HOW TO SWITCH BETWEEN BRUSING MODES:

- Power On/Off and brushing mode button(c).keep pressing this button,the operating sequence is power on,in Clean mode then White mode,then Polish mode,then Massage mode,then Sensitive mode ,end power off. The device will trun of if pressing any mode over 5 Sec.
- The RLT206 automatically starts in the Clean mode when pressing the brushing mode button once.
- The RLT206 starts from Clean mode,and switch to White,Polish,Massage and Sensitive by pressing on/off button different times.The device will turn off after Sensitive mode pressing.

HOW TO SWITCH BETWEEN BRUSING MODES:

- Attach a brush head to the handle, You can use any kind of toothpaste, To avoid splashing, guide the brush-head to your teeth before switching the appliance on.
- Place the bristles against the teeth at a slight angle towards the gum line Press the button to turn the toothbrush ,There are five cleaning modes .
- Take out an brush head (a) and the handle ,attach a brush head to the handle,there will be about 0.5 mm space between the brush head and the handle.
- Place the bristles against the teeth at a slight angle towards the gum line Press the button to turn the toothbrush .There are five cleaning modes.(Fig1)
- After brushing ,take down the brush head from the handle.(Fig2)
- Rinse both the brush head and handle under running water and wipe them dry.

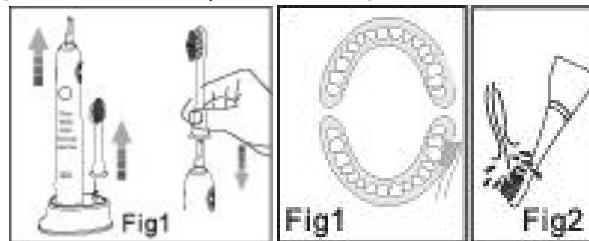

CLEANING RECOMMENDATIONS:

- After brushing, rinse your brush head thoroughly under running water.
- Remove the brush head from the handle and clean both parts separately under running water. Then wipe them dry before you put the handle onto the charging unit.
- Before cleaning your charging unit, make sure you unplug the adaptor. Only use a dry or damp-dry cloth to wipe the surface of the charging unit.

important: Do not clean the charging unit by immersing it under running or still water.

- Place the charging unit onto a dry and safe place.
- Do not plug in the charging unit until it is completely dry.

This product contains rechargeable batteries in the interest of protecting the environment.